



MDF-DC102VH-PE / MDF-DC202VH-PE

VIP ECO

-80°C Chest Freezers

84L / 180L



PHCbi's MDF-DC102VH-PE and MDF-DC202VH-PE VIP ECO Chest Freezers provide maximum sample storage capacity within an optimum footprint combined with industry leading cooling performance and reliability for an ideal long-term preservation solution.+

Optimum Sample Storage


The use of patented VIP PLUS panels within the freezer cabinet provides a reduced wall thickness for maximum interior volume. This technology allows the MDF-DC102VH-PE and MDF-DC202VH-PE to require a small installation space with a great storage capacity.

Efficient Refrigeration

Naturally occurring hydrocarbon (HC) refrigerants provide more efficient cooling due to their high latent heat of evaporation. As well as improved performance this leads to lower GWP.

Ease of Use & Intelligent Security

The OLED panel has good visibility and intuitive operation. The USB port and data log functions simplifies temperature data and events management

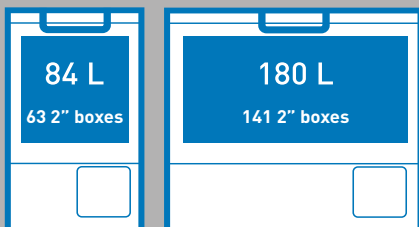

Personal Storage

Ideal for reliable storage of small volume clinical trial samples and temperature control testing as well as for a wide variety of other biological applications.



Dependable Preservation

Quality of design and manufacture ensures trusted and reliable storage that maintains the integrity of precious samples.



Uniform Sample Storage

The uniform and stable temperatures provided are ideal for storage of DNA, RNA and other important biological samples.

VIP ECO -80°C Chest Freezers



MDF-DC102VH-PE

MDF-DC202VH-PE

Specially developed single-compressor system

The unique auto-cascade system uses a single compressor for efficient cooling, enabling an extremely small footprint.

The highly efficient refrigeration system provides uniform and stable ultra-low temperatures down to -86°C.


Filterless Design

The unique cooling circuit makes the inconvenient customer maintenance procedure of filter cleaning unnecessary.

Minimal Footprint

PHCbi VIP ECO ULT freezers provide up to 30% more storage capacity than a conventionally insulated freezer without increasing the footprint.⁵⁾

The MDF-DC102VH-PE and MDF-DC202VH-PE feature VIP PLUS vacuum panels for increased energy efficiency and cooling performance.

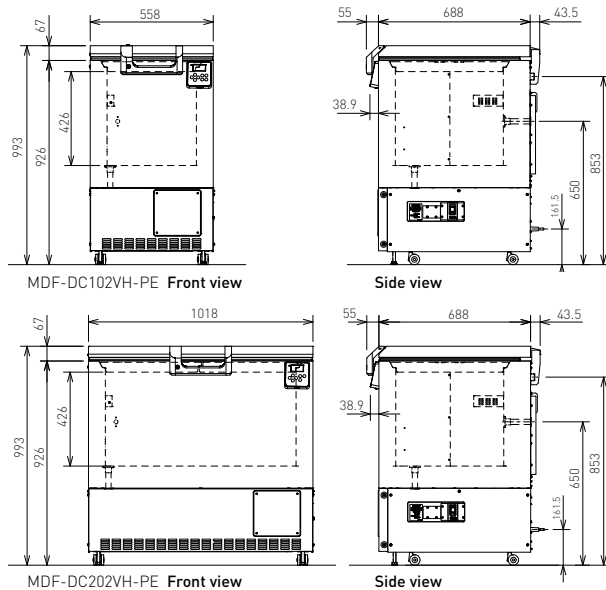

Model Number		MDF-DC102VH-PE	MDF-DC202VH-PE
Dimensions			
External dimensions (WxDxH) ¹⁾	mm	558 x 688 x 993	1018 x 688 x 993
Internal dimensions (WxDxH)	mm	405 x 490 x 426	865 x 490 x 426
Volume	litres	84	180
Net weight	kg	85	115
Capacity Vertical rack: NIR-207C	2" boxes	42	119
Capacity Ladder racks: VC-660/3SC	2" boxes	63	141
Performance			
Cooling performance ²⁾	°C	-86	
Temperature setting range	°C	-40 to -90	
Temperature control range ²⁾	°C	-40 to -86	
Control			
Controller	Microcomputer control system		
Display	OLED		
Temperature sensor	Pt1000		
Refrigeration			
Compressor	W	550	
Refrigerant	HC mixed refrigerant		
Refrigerant weight	g	80	110
Refrigerant GWP		5	6
Total Refrigerant weight (CO ₂ equivalent)	t	0.0004	0.0007
Insulation material	PUF / VIP Plus		
Insulation thickness	mm	76.5	
Construction			
Exterior material	Painted Steel		
Interior material	Painted Steel		
Outer Lid	qty	1	
Outer Lid Lock	Y		
Inner Lid	qty	1 (Styrofoam)	2 (Styrofoam)
Access Port	qty	2	
- position	Back/Bottom left front		
- diameter	∅ mm	17	
Casters	qty	4 (2 leveling feet)	
Data logging/USB	Y		
Alarms		(R = Remote Alarm, V = Visual Alarm, B = Buzzer Alarm)	
Power failure	V-B-R		
High temperature	V-B-R		
Low temperature	V-B-R		
Filter	Filterless design		
Electrical and Noise Level			
Power Supply	V	220/230/240	
Frequency	Hz	50	
Noise Level ³⁾	dB(A)	52	
Options			
Liquid CO ₂ back-up	MDF-UB8-PW		
Temperature recorders			
- Circular type ⁴⁾	MTR-G85C-PE		
- Chart paper	RP-G85-PW		
- Ink pen	PG-R-PW		
- Continuous strip type	MTR-85H-PW		
- Chart paper	RP-85-PW		
- Ink pen	DF-38FP-PW		
- Recorder housing	MDF-S3085-PW		

Appearance and specifications are subject to change without notice.

1) Exterior dimensions of main cabinet only, excluding handle and other external projections - See dimensions drawings on website for full details.

2) Air temperature measured at freezer centre, ambient temperature +30°C, no load.

3) Nominal value - Background noise 20dB(A).

4) Requires sensor cover MTR-C8-PW.

5) Compared with the MDF-193.

PHC Europe

A Member of PHC Group

Eikdonk 1 | 4825 AZ Breda | Netherlands
T: +31 (0) 76 543 3833

www.phcd.com/eu/biomedical